



RESMED INSTRUCTIONS

CHARGE THE BATTERY FIRST

- Connect the wall plug to the charger unit and fully charge the battery before first use.
- Plug one end of the charger into the wall and the other end to the battery pack.
- If the light in the corner of the rectangular box on charger is red, it is still charging. Once the light turns green, it is fully charged. Note: The red light on the cigarette lighter male end of the charger that connects to the battery will not turn green.



Female End of Battery

Wall Plug

Cigarette Lighter Male End

Light on the other side of the charger will turn green when fully charged

USING THE BATTERY WITH YOUR CPAP DEVICE

- When using the CPAP Battery with your CPAP device, use the appropriate DC cord that came with your battery.
- Plug the DC cord cigarette lighter male end to the battery and the other end of the DC cord into your machine.
- Turn on your machine and use as normal.

DC cord cigarette lighter male end



connect this end to the back of your ResMed CPAP machine



SAFETY & HANDLING OF BATTERY

- Battery should be charged prior to use.
- Use only the charger that came with the battery, do not try to charge it with an alternate device.
- Do not short circuit the battery, permanent damage to the battery could happen.
- Do not incinerate or mutilate the battery as it may cause it to burst or release toxic material.
- Avoid the battery being used in an airtight compartment. The battery should be ventilated.
- If you find any extremely high temperature or leakage, stop use.
- Allow battery to cool before charging.
- Do not remove the outer sleeve from the battery or cut into its housing.
- Disconnect the battery from your device when not in use.
- During long term storage the battery should be charged and discharged once every 6 months.
- Do not mix the battery with other battery brands.
- Never put the battery into water or seawater.
- Do not attempt to take the battery apart or puncture it. The alkaline electrolyte is harmful to eyes and skin.
- Keep away from children. If swallowed, contact a physician at once.

RESPIRONICS INSTRUCTIONS



CHARGE THE BATTERY FIRST

- Connect the wall plug to the charger unit and fully charge the battery before first use.
- Plug one end of the charger into the wall and the other end to the battery pack.
- If the light in the corner of the rectangular box on charger is red, it is still charging. Once the light turns green, it is fully charged. Note: The red light on the cigarette lighter male end of the charger that connects to the battery will not turn green.



Female End of Battery

Wall Plug


Cigarette Lighter Male End

Light on the other side of the charger will turn green when fully charged

USING THE BATTERY WITH YOUR CPAP DEVICE

- When using the CPAP Battery with your CPAP device, use the appropriate DC cord that came with your battery.
- Plug the DC cord cigarette lighter male end to the battery and the other end of the DC cord into your machine.
- Turn on your machine and use as normal.

DC cord cigarette lighter male end

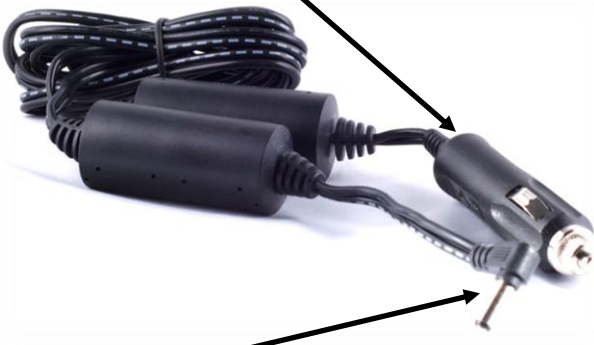

connect this end to the back of your Respiration CPAP machine



SAFETY & HANDLING OF BATTERY

- Battery should be charged prior to use.
- Use only the charger that came with the battery, do not try to charge it with an alternate device.
- Do not short circuit the battery, permanent damage to the battery could happen.
- Do not incinerate or mutilate the battery as it may cause it to burst or release toxic material.
- Avoid the battery being used in an airtight compartment. The battery should be ventilated.
- If you find any extremely high temperature or leakage, stop use.
- Allow battery to cool before charging.
- Do not remove the outer sleeve from the battery or cut into its housing.
- Disconnect the battery from your device when not in use.
- During long term storage the battery should be charged and discharged once every 6 months.
- Do not mix the battery with other battery brands.
- Never put the battery into water or seawater.
- Do not attempt to take the battery apart or puncture it, the alkaline electrolyte is harmful to eyes and skin.
- Keep away from children. If swallowed, contact a physician at once.